
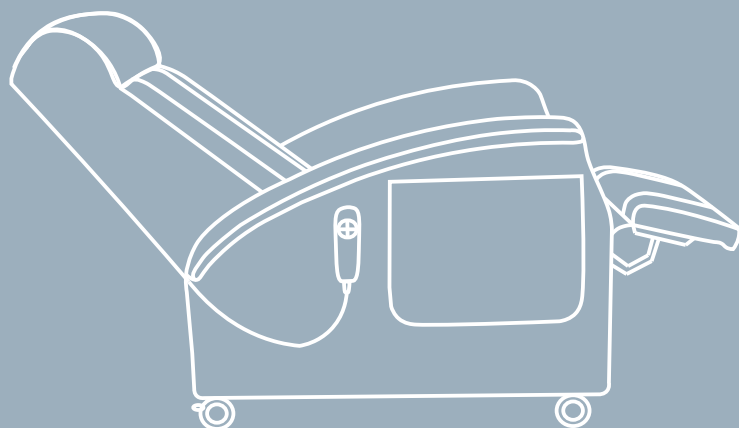


teal

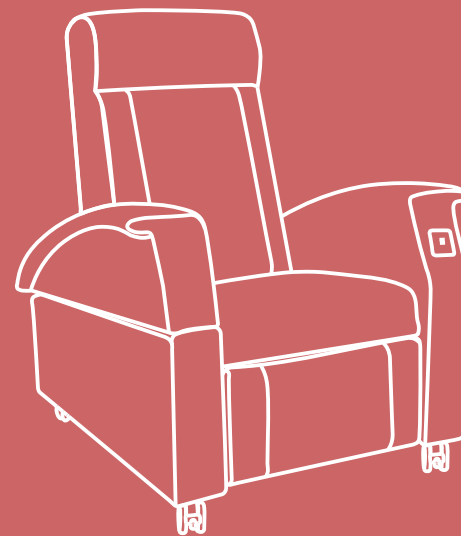

**Astra Stellar**  
Specialist Electric Recliner

# Astra Stellar specialist electric recliner.


Designed with enhanced patient features, for treatment areas including Oncology, Renal and Recovery Suites.





**Astra Stellar has a high specification as standard including emergency back release and can be customised from a suite of options.**



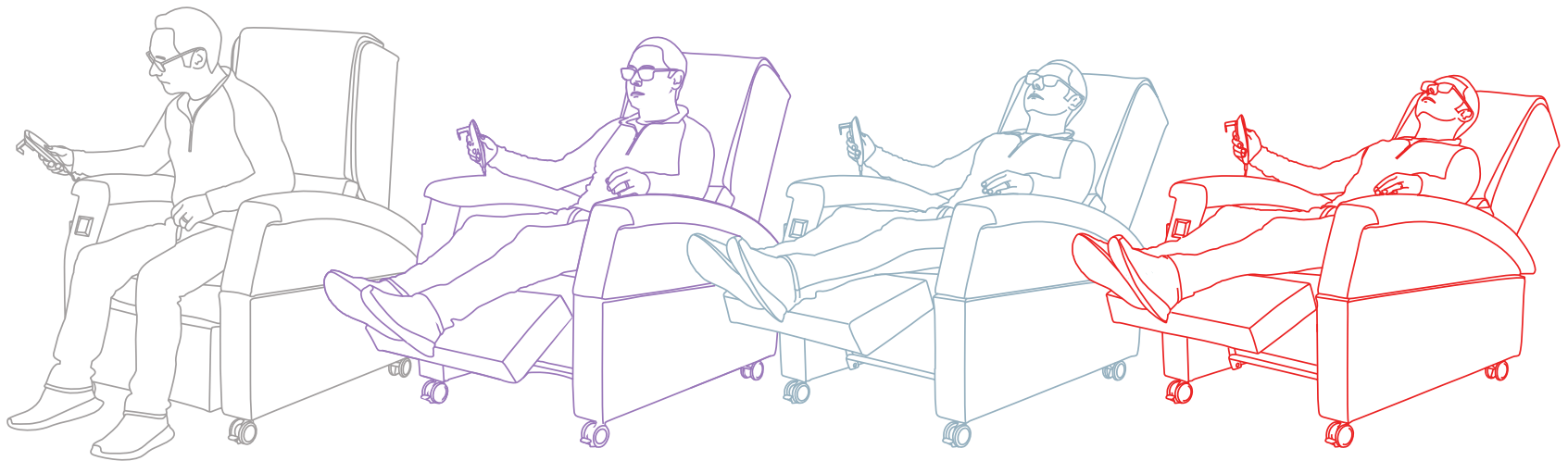

**Handset operated tilt-in-space electric recline with dual motor allowing independent back and footrest operation.**

Rise on REST22 model.




# Features

## Standard Features

- Postural support with Reflexion® seat, arm tops and leg support as standard for comfort and pressure management
- Tilt-in-space functionality as standard (backrest and seat recline synchronously, helping to maintain posture)
- Extra wide 150mm arms for added comfort during procedures
- Emergency quick release back
- Dual motor allows independent back and footrest operation
- Recline angle 128° (REST21 Recline)
- Recline angle 158° (REST22 Rise and Recline)
- Manoeuvring handle for ease of movement (not for portering)
- Four heavy duty 65mm lockable castors
- Max load: 125kg (20st) REST21
- Max load: 160kg (25st) REST22

## Optional Features

- Removable arm(s) for patient transfer
- Waterfall back
- Waterfall back with containment wings
- Anti-pinch leg support
- Lithium Ion battery pack (100 operations)
- USB charging for Iphone/Ipad
- Folding table (anti-microbial polymer)

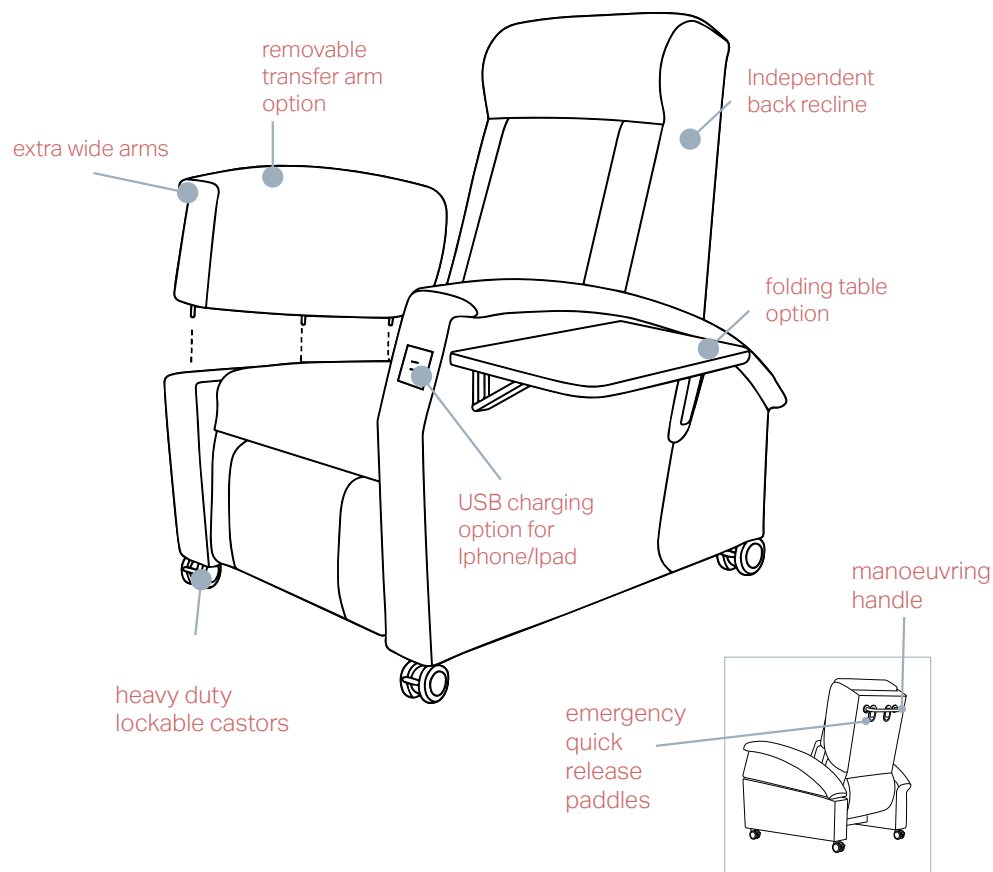

## Model Options

**REST21 Electric Recliner**

**REST22 Electric Riser/Recliner**



## Emergency quick release back



1. Ensure the rear braked castors are in the locked position

2. With both hands on the manoeuvring handle on the rear of the back rest, place two fingers under each paddle.



3. Push the back slightly forwards while pulling the paddles towards you. This will allow the back to release and lower quickly.

4. To return to an upright position push back up from both corners.



# Optional features



1. standard back 2. waterfall back with wings 3. removable transfer arm 4. footrest stitching (with standard back)  
5. folding table 6. usb charger



# Technical Specifications

## Model Options:

### REST21 Electric Recliner


Recline angle: 128°  
Max load 125kg (20st)

### REST22 Electric Riser/Recliner

Recline angle: 158°  
Max load 160kg (25st)




BB pressure relief foam

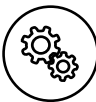

Dual Motor



Tilt in Space

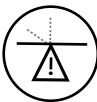

Emergency Back Relief



Braked Castors



Removable Arm Option




Twin USB Port Option




Lithium Battery Pack Option




Folding Table Option

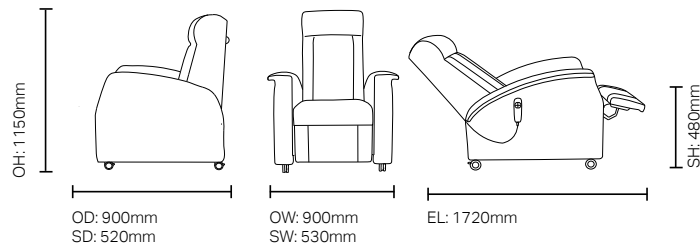

## Both Models:

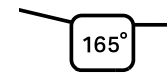

Emergency quick release back  
Recline Angle 165°

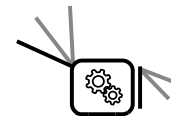

Dual Motor allows for independent  
backrest and footrest recline

## Environmental Product Analysis

Recyclability: 99%  
Recycled Content: 45.15%





# Clatterbridge Hospital Case Study

Teal HealthCare helps deliver world class cancer facilities at the new Liverpool Clatterbridge Hospital

A spectacular new hospital delivering highly-specialist care including pioneering immunotherapy and the most advanced forms of radiotherapy has been opened in the heart of Liverpool, providing world-class cancer facilities for the community in Cheshire and Merseyside.



Astra Stellar in use at Clatterbridge Hospital

Clatterbridge Cancer Centre - Liverpool is part of a £162 million project to transform cancer services across the Cheshire and Merseyside region. Working alongside key partners including clinical staff, NHS Shared Business Services and the main contractor Laing O'Rourke, Teal HealthCare installed Astra Stellar recliners, Astra Club recliners, patient seating and specialist healthcare products throughout the spectacular new 11-storey facility.



# teal

Teal  
Branch Road  
Lower Darwen  
Lancashire  
BB3 0PR

Tel: +44 (0) 1254 688210  
Sales Freephone: 0800 980 8998  
Email: [sales@teal.co.uk](mailto:sales@teal.co.uk)  
[teal.co.uk](http://teal.co.uk)

**The Senator Group**

## Contact

**Lisa Strong**  
Area Business Manager - South East  
[lisa@teal.co.uk](mailto:lisa@teal.co.uk)  
+44 (0) 7718 660401

**Jane Staufenbiel**  
Area Business Manager - Midlands  
[jane@teal.co.uk](mailto:jane@teal.co.uk)  
+44 (0) 7725 273149

**Ben Wooster**  
Area Business Manager -  
South West & Wales  
[ben@teal.co.uk](mailto:ben@teal.co.uk)  
+44 (0) 7725 273148

**Carmel Satterthwaite**  
Area Business Manager - North of England,  
North Wales, N.I & Ireland  
[carmel@teal.co.uk](mailto:carmel@teal.co.uk)  
+44 (0) 7501 254330

**Ciaran Claxton**  
Business Development Manager -  
R.O.I & Northern Ireland  
[Ciaran@teal.co.uk](mailto:Ciaran@teal.co.uk)  
+44 (0) 7836 366 536

**Scotland - Langstane Press**  
[www.langstane.co.uk](http://www.langstane.co.uk)  
[John.Lister@langstane.co.uk](mailto:John.Lister@langstane.co.uk)  
01506 674520

