

**PSTL01, PSTL02, PSTL03, PSTL04, PSTL05,  
PSTL06, PSTL07, PSTL08**

The Tolero metal patient chair demands attention by virtue of its slimline profile, contemporary form and a host of unrivalled features for infection control, tissue viability and manual handling, including: removable seat, Reflexion® high risk pressure distributing seat, housekeeping wheels and integral height adjustable legs (**PSTL02 / 03 / 05 / 06 models only**).



### Standard Features


- Strong, durable high grade steel frame creates a modern aesthetic
- Modern, slim profile with contoured back mimics curvature of spine for constant support
- Elegant moulded black plastic arms and extruded frame with integral anti-bacterial technology
- Removable depth adjustable seat (2 positions) for cleaning / access
- Easy clean polymer moulded seat base
- Reflexion® BB medium / high risk pressure management seat
- Sacral gap as standard for cleaning/ infection management
- Ergonomically positioned arms for easy sit to stand transition
- Integral housekeeping wheels as standard
- Height adjustable legs up to 100mm (25mm increments) (PSTL02 / 03 / 05 / 06 models only)
- Resilient to commonly used healthcare cleaning systems
- Maximum user weight: 160kg (25st)  
Bariatric user weight 320kg (50st)
- Black arms as standard (All models)

### Optional Features


- Frame finish options (Stone, Onyx, Anthracite)
- Fixed seat

## Technical Information


### Dimensions:

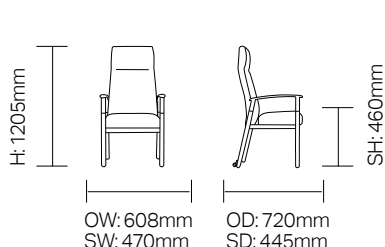

**PSTL01**  
High Back Chair  
with Housekeeping Wheels

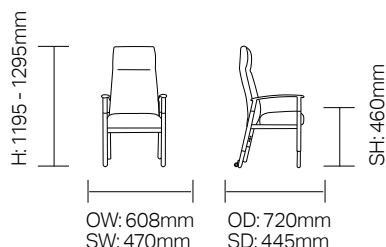


**PSTL02 Tolero**  
High Back Chair,  
Integral Height Adjustable Legs  
with Housekeeping Wheels

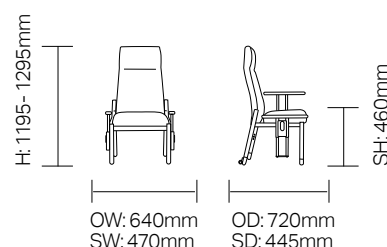


**PSTL03 Tolero**  
High Back Chair, Drop Arms,  
Integral Height Adjustable Legs  
with Housekeeping Wheels


## Dimensions:

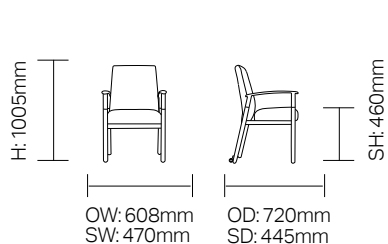


**PSTL04**  
Tolero Mid Back Chair  
with Housekeeping Wheels

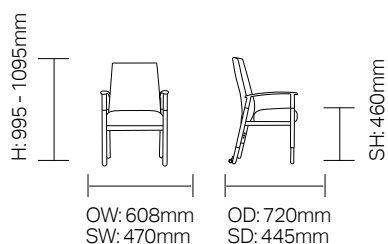


**PSTL05**  
Tolero Mid Back Chair,  
Height Adjustable Legs  
with Housekeeping Wheels

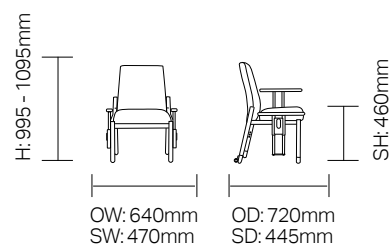

**PSTL06**  
Tolero Mid Back Chair, Drop Arms,  
Height Adjustable Legs  
with Housekeeping Wheels

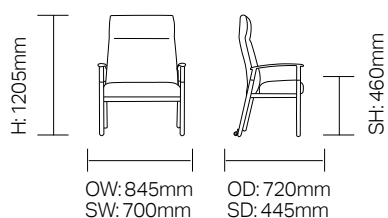

**PSTL07 Tolero**  
High Back Bariatric Chair  
with Housekeeping Wheels

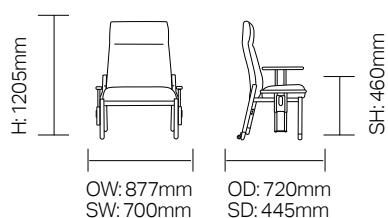

**PSTL08 Tolero**  
High Back Bariatric Chair,  
Drop Arms with Housekeeping  
Wheels



## Tolero Features




### Removable Seat

Removable seat with polymer base facilitates cleaning and infection control.



### Adjustable Seat Depth

Adjustable seat allows 2 depths ensuring support from the chair back without cutting off circulation at the knees.

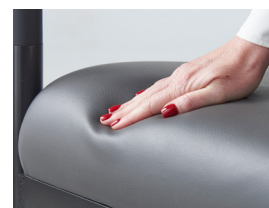

### Pressure Care Seat

Reflexion® BB medium / high risk pressure management seat. Standard on all models.

Reflexion® CC high risk pressure management seat option also available.



### Height Adjustable Legs (Model Specific)

A range of chair heights provide appropriate seating height and bring benefit to posture, circulation, balance on rising and core strength. Image shows maximum height. In 25mm increments.



### Ergonomic slim profiled back

Slim profile with contoured moulded back mimics curvature of spine for constant support.



### Easy Rise Arms

The ability to stand from a sitting position is key to maintaining mobility and independence and is an important factor for patients.



### Sacral Gap

Sacral gap for ease of cleaning & access.



### Housekeeping Wheels

Housekeeping wheels facilitate manual handling. Housekeeping wheels standard on all models.

**Strong, durable high grade tubular metal frame with anti-bacterial easy clean coated finish for infection management.**

## Finish Options

Frame colour options:



Stone



Onyx



Anthracite

Colour Disclaimer: Printed images and digital images may appear differently to the final finish.

Upholstery Selections can be found [Here](#)



**PSTL08** Tolero Bariatric Patient Chair shown above



### 3 R's

Teal is committed to continually improving the sustainability of all environmental aspects within our business.

To meet both international standards and our own environmental targets we apply the three R's principle

## Reduce, Reuse and Recycle.

Whilst recycling is the element which receives the most exposure it is actually the last option available and should never be the prime target in anyone's battle to reduce waste.

It is our duty as individuals and as a company to initially attempt to Reduce usage. Then we should look to Reuse wherever possible and finally, only after these two processes have been exhausted, should we consider Recycling.



**Life cycle understanding**



**Waste elimination**



**Recycled content innovation**



**End of life solutions**

teal

Branch Road  
Lower Darwen  
Lancashire UK  
BB3 0PR

Tel: +44 (0) 1254 688210  
Sales Freephone: 0800 980 8998  
Email: [sales@teal.co.uk](mailto:sales@teal.co.uk)  
[www.teal.co.uk](http://www.teal.co.uk)

The Senator Group

